

3-アミノフェノールのマウスを用いた  
経口投与によるがん原性試験（混水試験）報告書

試験番号：0712

# APPENDICES

## APPENDICES

- APPENDIX 1-1 IDENTITY OF 3-AMINOPHENOL IN THE 2-YEAR DRINKING WATER STUDY
- APPENDIX 1-2 STABILITY OF 3-AMINOPHENOL IN THE 2-YEAR DRINKING WATER STUDY
- APPENDIX 1-3 CONCENTRATION OF 3-AMINOPHENOL IN FORMULATED WATER IN THE 2-YEAR DRINKING WATER STUDY
- APPENDIX 1-4 STABILITY OF 3-AMINOPHENOL IN FORMULATED WATER
- APPENDIX 2 METHODS, UNITS AND DECIMAL PLACE FOR HEMATOLOGY AND BIOCHEMISTRY IN THE 2-YEAR DRINKING WATER STUDY OF 3-AMINOPHENOL

APPENDIX 1-1

IDENTITY OF 3-AMINOPHENOL  
IN THE 2-YEAR DRINKING WATER STUDY

# IDENTITY OF 3-AMINOPHENOL IN THE 2-YEAR DRINKING WATER STUDY

Test Substance : 3-Aminophenol (Wako Pure Chemical Industries, Ltd.)

A. Lot No. : WKE2284



## 1. Spectral Data

### Mass Spectrometry

Instrument : Agilent Technologies 5973N Mass Spectrometer

Ionization : EI (Electron Ionization)

Ionization Voltage : 70eV

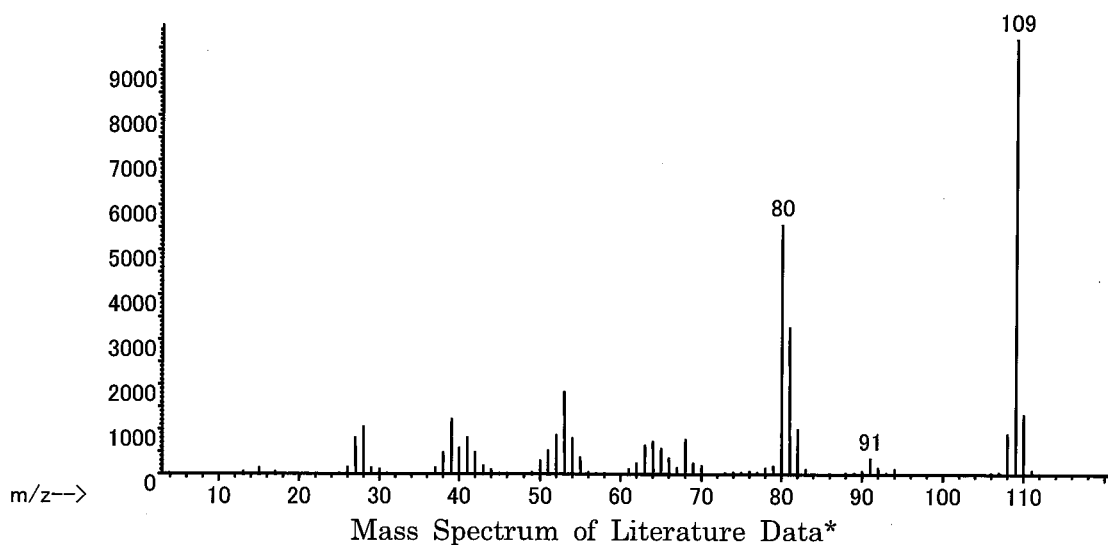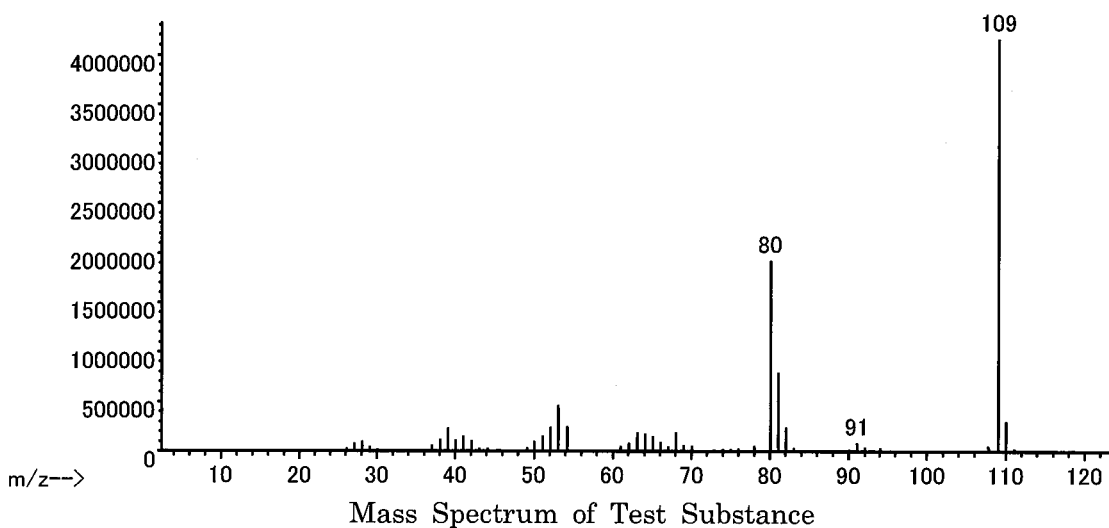


Result: The mass spectrum was consistent with literature spectrum.  
(\*McLafferty FW, ed. 1994. Wiley Registry of Mass Spectral Data. 6th ed.  
New York, NY : John Wiley and Sons.)

## Infrared Spectrometry


Instrument : Shimadzu FTIR-8200PC Infrared Spectrometer

Cell : KBr

Resolution : 4  $\text{cm}^{-1}$

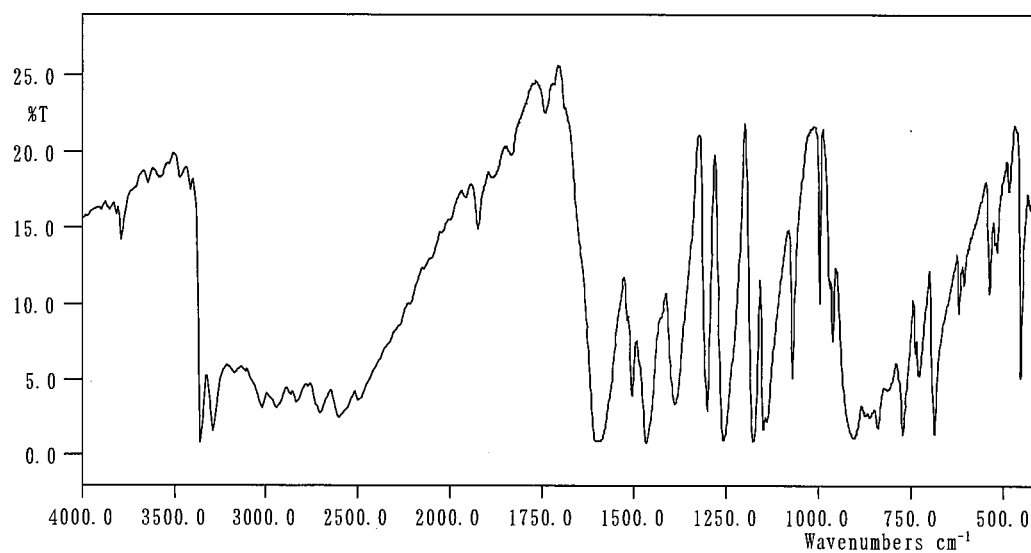

Infrared Spectrum of Test Substance

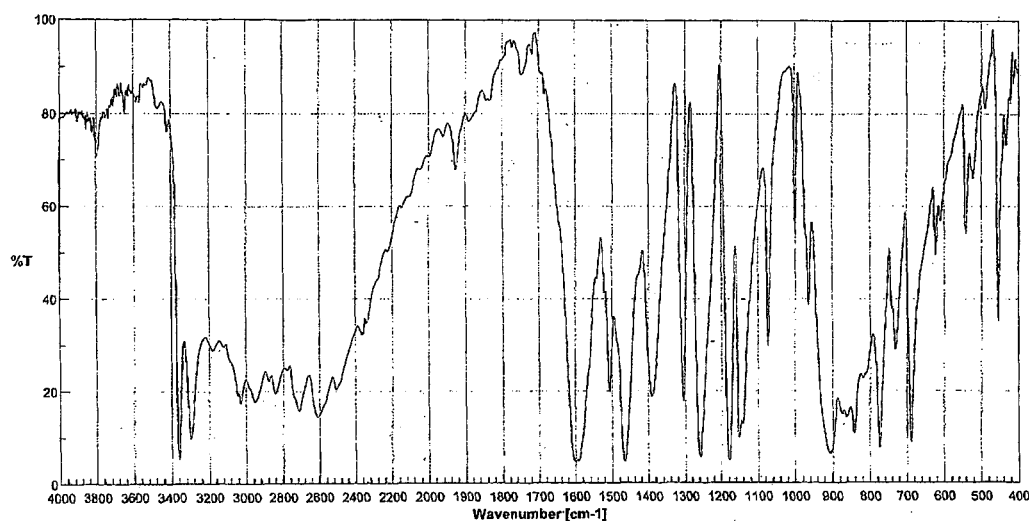

Infrared Spectrum of Literature Data\*

Result: The infrared spectrum was consistent with literature spectrum.  
(\*Performed by Wako Pure Chemical Industries, Ltd.)

2. Conclusion: The test substance was identified as 3-aminophenol by mass spectrum and infrared spectrum.

B. Lot No. : CDQ4568



1. Spectral Data

Mass Spectrometry

Instrument : Agilent Technologies 5973N Mass Spectrometer

Ionization : EI (Electron Ionization)

Ionization Voltage : 70eV

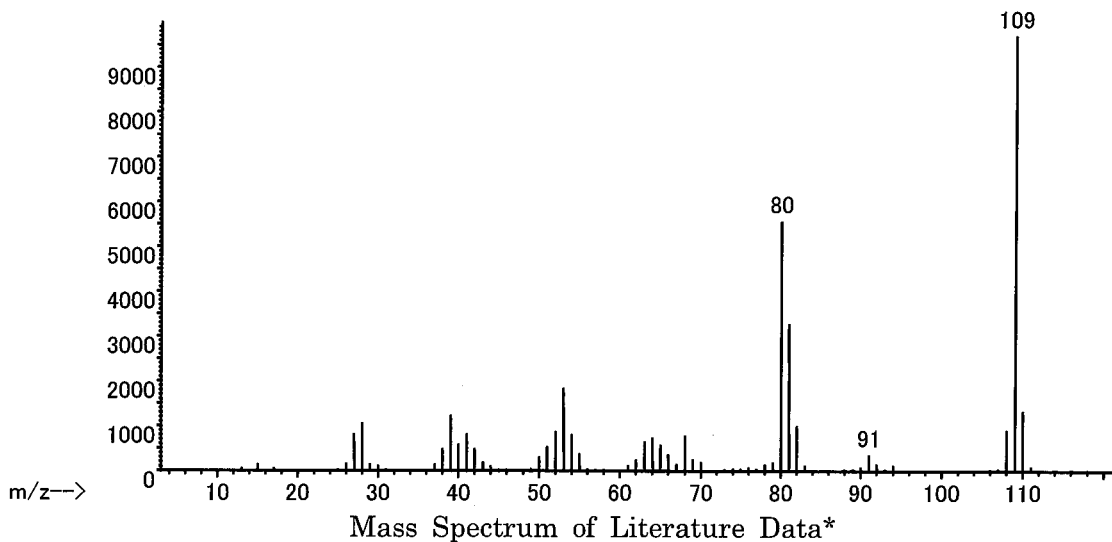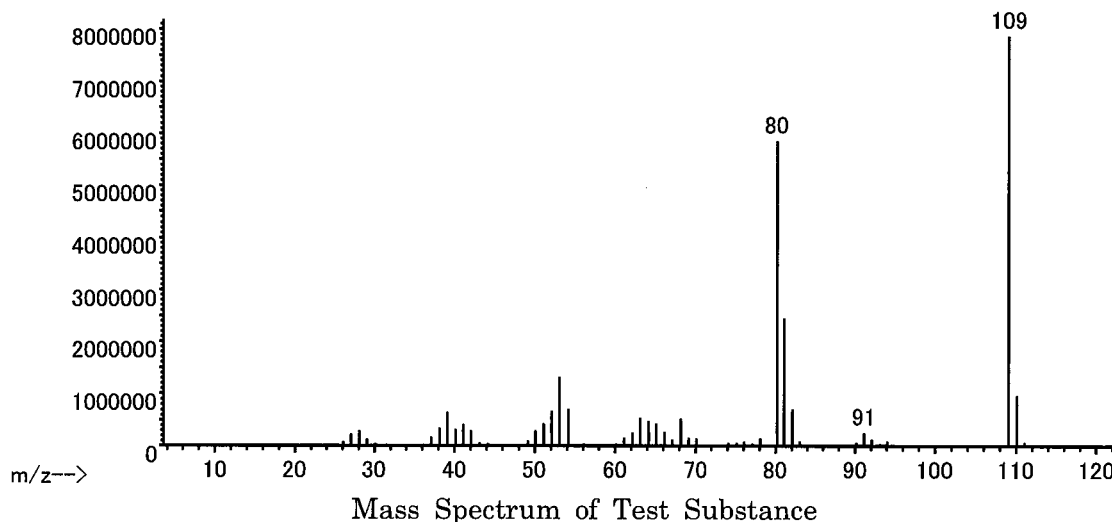


Result: The mass spectrum was consistent with literature spectrum.  
(\*McLafferty FW, ed. 1994. Wiley Registry of Mass Spectral Data. 6th ed.  
New York, NY : John Wiley and Sons.)

## Infrared Spectrometry


Instrument : Shimadzu FTIR-8200PC Infrared Spectrometer

Cell : KBr

Resolution : 4  $\text{cm}^{-1}$

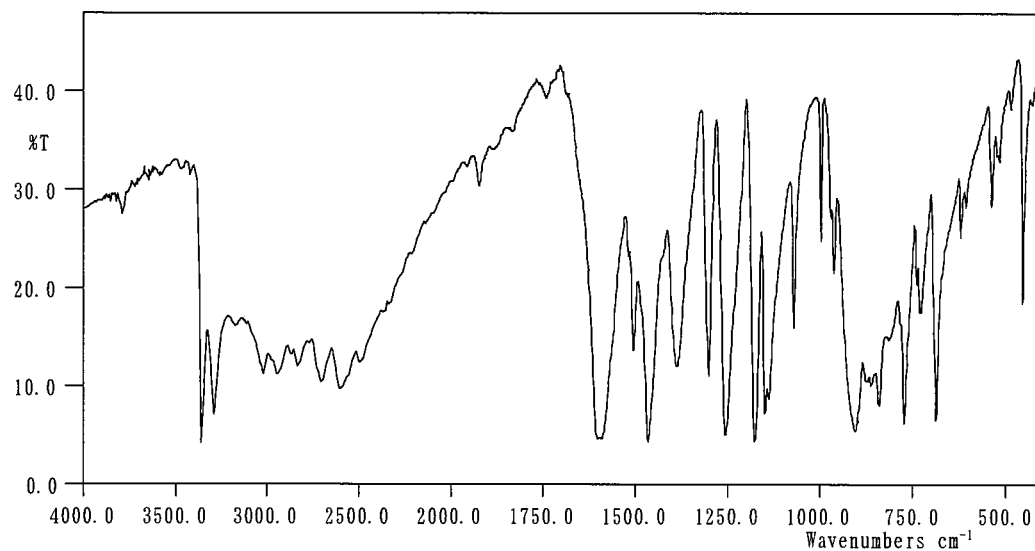

Infrared Spectrum of Test Substance

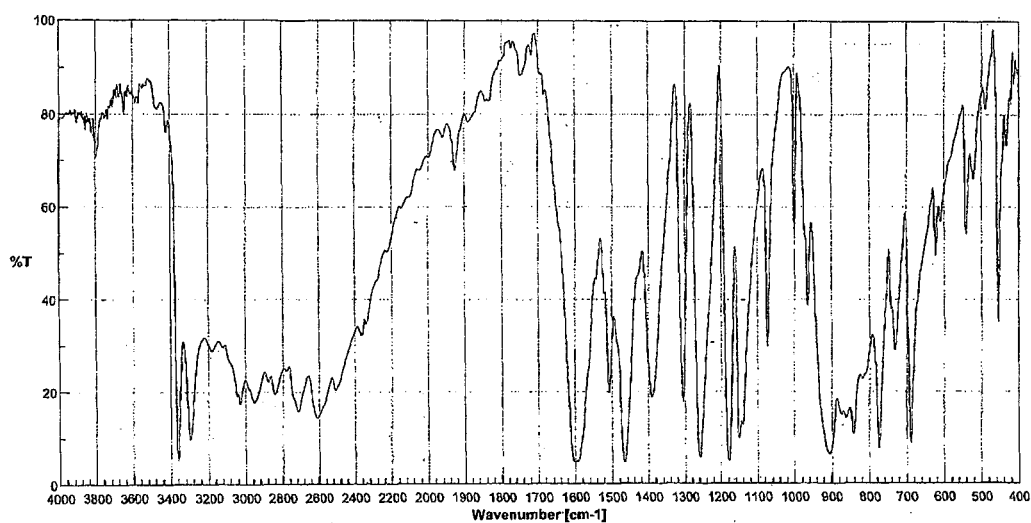

Infrared Spectrum of Literature Data\*

Result: The infrared spectrum was consistent with literature spectrum.  
(\*Performed by Wako Pure Chemical Industries, Ltd.)

2. Conclusion: The test substance was identified as 3-aminophenol by mass spectrum and infrared spectrum.

APPENDIX 1-2

STABILITY OF 3-AMINOPHENOL  
IN THE 2-YEAR DRINKING WATER STUDY



## STABILITY OF 3-AMINOPHENOL IN THE 2-YEAR DRINKING WATER STUDY

Test Substance : 3-Aminophenol (Wako Pure Chemical Industries, Ltd.)

A. Lot No. : WKE2284

## 1. High Performance Liquid Chromatography

Instrument : Shimadzu LC-10 High Performance Liquid Chromatograph

Column : TSK-GEL ODS-80TM (4.6 mm  $\phi$   $\times$  15 cm)

Column Temperature: 40 °C

Flow Rate : 1 mL/min

Mobile Phase : Acetonitrile : Methanol : 5mM Sodium dodecyl sulfate solution (phosphoric acid pH2.2)  
= 3 : 3 : 4

Detector : UV (275 nm)

Injection Volume : 10  $\mu$ L

Date Analyzed	Peak No.	Retention Time (min)	Area (%)
2008.07.02	1	3.960	100
2009.11.19	1	3.949	100

Result: High performance liquid chromatography indicated one major peak (peak No.1) analyzed on 2008.7.2 and one major peak (peak No.1) analyzed on 2009.11.19. No new trace impurity peak in the test substance analyzed on 2009.11.19 was detected.

2. Conclusion: The test substance was stable for the period that the test substance had been used for the study.

B. Lot No. : CDQ4568

1. High Performance Liquid Chromatography

Instrument : Shimadzu LC-10 High Performance Liquid Chromatograph

Column : TSK-GEL ODS-80TM (4.6 mm  $\phi$   $\times$  15 cm)

Column Temperature: 40 °C

Flow Rate : 1 mL/min

Mobile Phase : Acetonitrile : Methanol : 5mM Sodium dodecyl sulfate solution (phosphoric acid pH2.2)  
= 3 : 3 : 4

Detector : UV (275 nm)

Injection Volume : 10  $\mu$ L

Date Analyzed	Peak No.	Retention Time (min)	Area (%)
2009.11.06	1	3.941	100
2010.08.02	1	3.967	100

Result: High performance liquid chromatography indicated one major peak (peak No.1) analyzed on 2009.11.6 and one major peak (peak No.1) analyzed on 2010.8.2. No new trace impurity peak in the test substance analyzed on 2010.8.2 was detected.

2. Conclusion: The test substance was stable for the period that the test substance had been used for the study.

APPENDIX 1-3

CONCENTRATION OF 3-AMINOPHENOL  
IN FORMULATED WATER IN THE 2-YEAR DRINKING WATER  
STUDY

CONCENTRATION OF 3-AMINOPHENOL IN FORMULATED WATER IN THE 2-YEAR DRINKING WATER STUDY

Analytical Method : The samples were analyzed by high performance liquid chromatography.

Instrument : Shimadzu LC-10 High Performance Liquid Chromatograph

Column : TSK-GEL ODS-80TM (4.6 mm  $\phi$   $\times$  15 cm)

Column Temperature: 40 °C

Flow Rate : 1 mL/min

Mobile Phase : Acetonitrile : Methanol : 5mM Sodium dodecyl sulfate solution (phosphoric acid pH2.2)  
= 3 : 3 : 4

Detector : UV (275 nm)

Injection Volume : 10  $\mu$ L

Date Analyzed	Target Concentration		
	625 <sup>a</sup>	1250	2500
2008.07.25	620 <sup>b</sup> (99.2) <sup>c</sup>	1250 (100)	2500 (100)
2008.09.26	637 <sup>b</sup> (102) <sup>c</sup>	1270 (102)	2560 (102)
2008.12.19	610 <sup>b</sup> (97.6) <sup>c</sup>	1320 (106)	2650 (106)
2009.03.13	592 <sup>b</sup> (94.7) <sup>c</sup>	1200 (96.0)	2470 (98.8)
2009.06.05	628 <sup>b</sup> (100) <sup>c</sup>	1260 (101)	2530 (101)
2009.08.28	613 <sup>b</sup> (98.1) <sup>c</sup>	1270 (102)	2550 (102)
2009.11.20	620 <sup>b</sup> (99.2) <sup>c</sup>	1220 (97.6)	2520 (101)
2010.02.12	640 <sup>b</sup> (102) <sup>c</sup>	1280 (102)	2550 (102)
2010.05.07	625 <sup>b</sup> (100) <sup>c</sup>	1250 (100)	2490 (99.6)

<sup>a</sup> ppm

<sup>b</sup> ppm (Mean measured concentration.)

<sup>c</sup> % (Mean measured concentration/target concentration  $\times$  100.)

APPENDIX 1-4

STABILITY OF 3-AMINOPHENOL IN FORMULATED WATER

## STABILITY OF 3-AMINOPHENOL IN FORMULATED WATER

Analytical Method : The samples were analyzed by high performance liquid chromatography.

Instrument : Shimadzu LC-10 High Performance Liquid Chromatograph

Column : TSK-GEL ODS-80TM (4.6 mm  $\phi$   $\times$  15 cm)

Column Temperature: 40 °C

Flow Rate : 1 mL/min

Mobile Phase : Acetonitrile : Methanol : 5mM Sodium dodecyl sulfate solution (phosphoric acid pH2.2)  
= 3 : 3 : 4

Detector : UV (275 nm)

Injection Volume : 10  $\mu$ L

Date Analyzed	Target Concentration	
	100 <sup>a</sup>	7500
2007.07.27	101 (100) <sup>b</sup>	7720 (100)
2007.07.31 <sup>c</sup>	97.5 ( 96.5)	7800 (101)

<sup>a</sup> ppm

<sup>b</sup> % (Percentage was based on the concentration on date of preparation.)

<sup>c</sup> Animal room samples

## APPENDIX 2

METHODS, UNITS AND DECIMAL PLACE FOR HEMATOLOGY  
AND BIOCHEMISTRY IN THE 2-YEAR DRINKING WATER  
STUDY OF 3-AMINOPHENOL

**METHODS, UNITS AND DECIMAL PLACE FOR HEMATOLOGY AND BIOCHEMISTRY  
IN THE 2-YEAR DRINKING WATER STUDY OF 3-AMINOPHENOL**

Item	Method	Unit	Decimal place
<b>Hematology</b>			
Red blood cell (RBC)	Light scattering method <sup>1)</sup>	$\times 10^9/\mu\text{L}$	2
Hemoglobin(Hgb)	Cyanmethemoglobin method <sup>1)</sup>	g/dL	1
Hematocrit(Hct)	Calculated as $\text{RBC} \times \text{MCV}/10$ <sup>1)</sup>	%	1
Mean corpuscular volume(MCV)	Light scattering method <sup>1)</sup>	fL	1
Mean corpuscular hemoglobin(MCH)	Calculated as $\text{Hgb}/\text{RBC} \times 10$ <sup>1)</sup>	pg	1
Mean corpuscular hemoglobin concentration (MCHC)	Calculated as $\text{Hgb}/\text{Hct} \times 100$ <sup>1)</sup>	g/dL	1
Platelet	Light scattering method <sup>1)</sup>	$\times 10^3/\mu\text{L}$	0
Reticulocyte	Light scattering method <sup>1)</sup>	%	1
Methemoglobin	Van Assendelft method <sup>2)</sup>	%	1
White blood cell(WBC)	Light scattering method <sup>1)</sup>	$\times 10^3/\mu\text{L}$	2
Differential WBC	Visual observation method (Wright staining)	%	0
<b>Biochemistry</b>			
Total protein(TP)	Biuret method <sup>3)</sup>	g/dL	1
Albumin (Alb)	BCG method <sup>3)</sup>	g/dL	1
A/G ratio	Calculated as $\text{Alb}/(\text{TP} - \text{Alb})$ <sup>3)</sup>	—	1
T-bilirubin	Azobilirubin method <sup>3)</sup>	mg/dL	2
Glucose	GlcK·G-6-PDH method <sup>3)</sup>	mg/dL	0
T-cholesterol	CE·COD·POD method <sup>3)</sup>	mg/dL	0
Triglyceride	MGLP·GK·GPO·POD method <sup>3)</sup>	mg/dL	0
Phospholipid	PLD·ChOD·POD method <sup>3)</sup>	mg/dL	0
Aspartate aminotransferase (AST)	JSCC method <sup>3)</sup>	U/L	0
Alanine aminotransferase (ALT)	JSCC method <sup>3)</sup>	U/L	0
Lactate dehydrogenase (LDH)	JSCC method <sup>3)</sup>	U/L	0
Alkaline phosphatase (ALP)	JSCC method <sup>3)</sup>	U/L	0
$\gamma$ -Glutamyl transpeptidase ( $\gamma$ -GTP)	JSCC method <sup>3)</sup>	U/L	0
Creatine kinase (CK)	JSCC method <sup>3)</sup>	U/L	0
Urea nitrogen	Urease·GLDH method <sup>3)</sup>	mg/dL	1
Sodium	Ion selective electrode method <sup>3)</sup>	mEq/L	0
Potassium	Ion selective electrode method <sup>3)</sup>	mEq/L	1
Chloride	Ion selective electrode method <sup>3)</sup>	mEq/L	0
Calcium	OCPC method <sup>3)</sup>	mg/dL	1
Inorganic phosphorus	PNP·XOD·POD method <sup>3)</sup>	mg/dL	1

1) Automatic blood cell analyzer (ADVIA120 : Siemens Healthcare Diagnostics Inc.)

2) Spectrophotometer (DU-530 : Beckman Coulter, Inc.)

3) Automatic analyzer (Hitachi 7080 : Hitachi, Ltd.)